





Complete Workstation




eXACT™ Temperature Control and Thermal Cycling

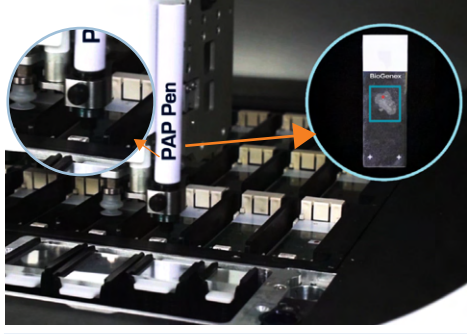


Hydrophobic Barrier

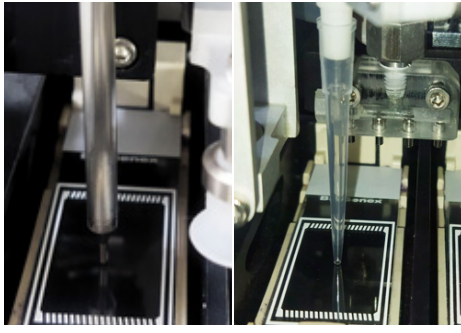


Reagent Dispensing

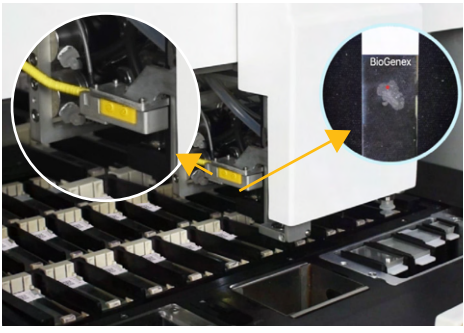

Reagent Dispensing at the designated spot

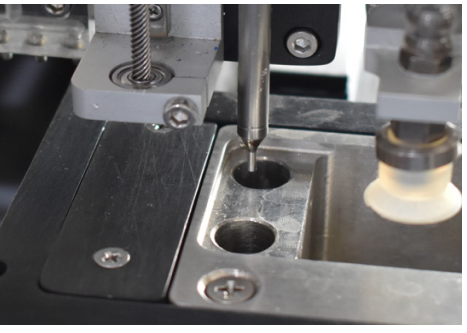

Probe Washing

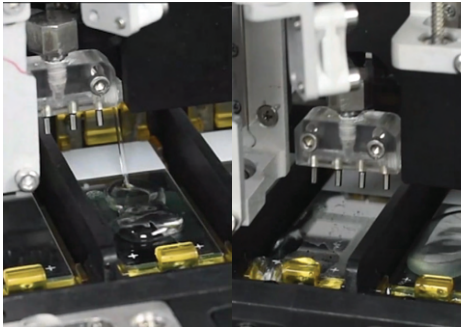


Wash & Air blow



Apply Coverslip

NanoVip™300

All-in-One FISH · ISH · IHC

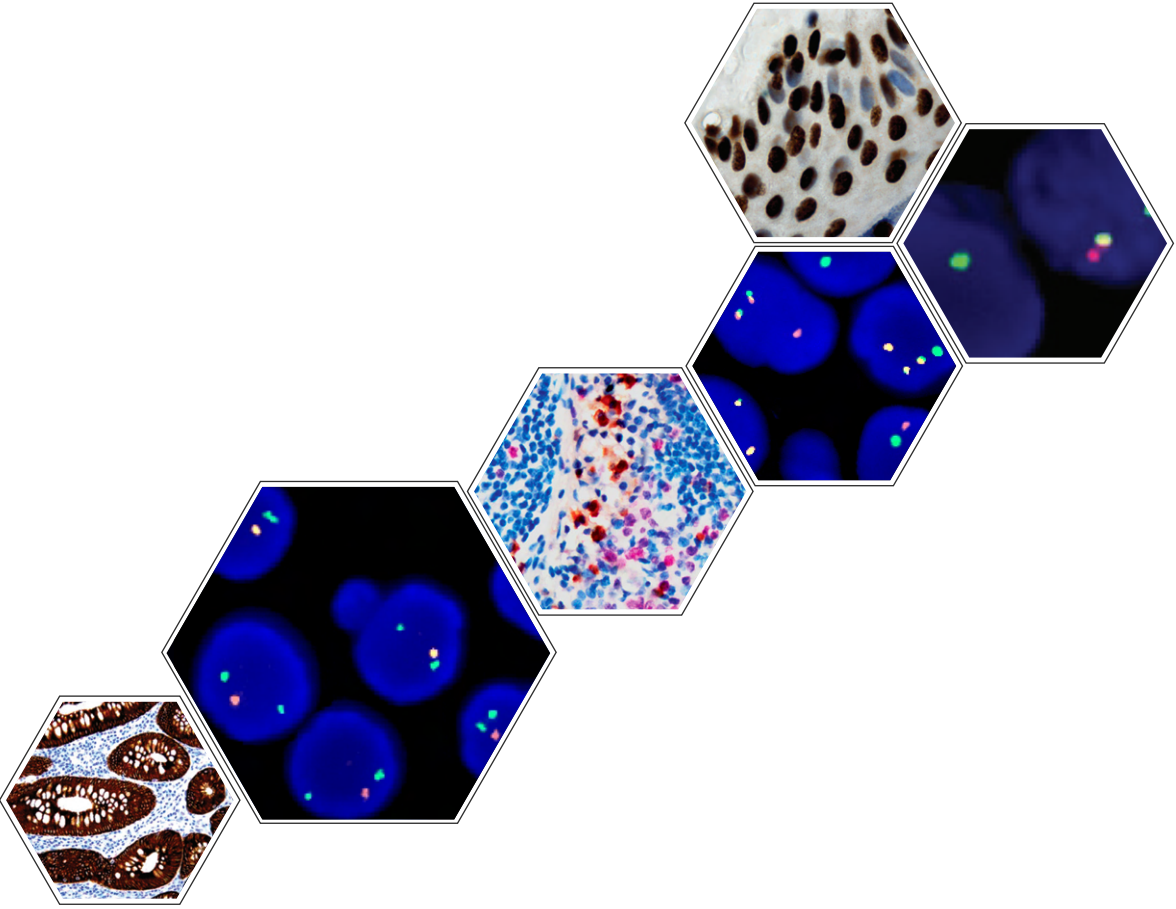

- Fully Automated Microtome to Microscope
- Automates any slide-based assay
- Baking to DAPI (FISH) & Final Coverslip
- In-built Camera detects the tissue size & location
- Draws hydrophobic barrier around the tissue section
- Dispenses as little as 2 µL at the selected spot
- Liquid level sensors for accurate liquid handling



Three Simple Steps:



Contact Information:

USA : 48810 Kato Road, Suite 100E & 200E, Fremont, CA 94538 USA
India : Aerospace & Precision Engineering, TSIIIC SEZ, Adibatla, Telangana-501510
China: Shanghai Exhibition Center, No.1000 Middle Yan An Road, Shang Hai Shi, China, 200041

Customer service - US : customerservice@biogenex.com
Customer service - India : indiacs@biogenex.com
Customer service - Global : internationalcs@biogenex.com
Instrument support : instrument.care@biogenex.com
Reagent support : reagentsupport@biogenex.com

For US	For India	For China	For Other Countries
T: +1-510-824-1400	T: +91-6309448856	T: +86-21-2220-6680	T: +91-7702301043
T: +1-800-421-4149		F: +86-21-2220-6679	
F: +1-510-824-1490			




www.biogenex.com


Video Link: <https://biogenex.com/product/nanovip-300/>

NanoVip™300


All-in-One FISH · ISH · IHC

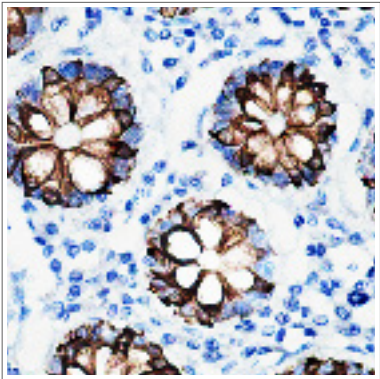

Immunohistochemistry

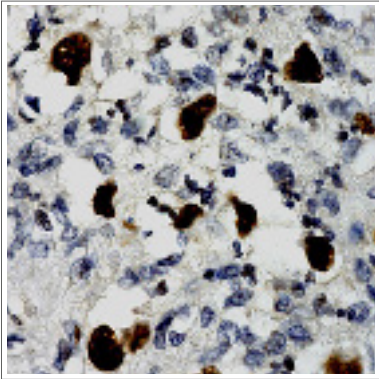

In Situ Hybridization

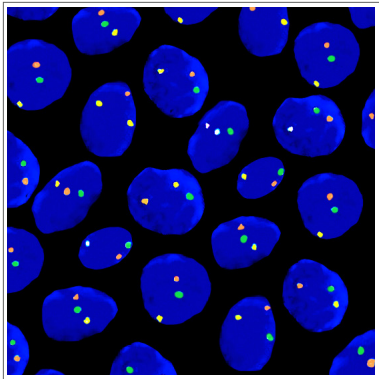

Fluorescence
In Situ Hybridization



NanoVip™ 300 System Specifications

Automation	Full (Baking through Final Coverslip)
Slide capacity	30 Standard Slides or Microchamber
Reagent Vials	2 x 20 Reagent Vials
Reagent vial capacity	22 mL
Bulk reagent	5 x 2 L
Waste container	1 x 10 L
Slide sensing	Open
Reagent sensing	Open
Reagent Dispense volumes	2 µL - 180 µL
Reagent Dispenser	1000 µL & 200 µL Pipette tips & Probe
Slide Temperature range	25 - 105 °C
Configuration	Bench Top
Protocols	Open System program any manual protocol
Operating Temperature	15 - 30 °C
Operating Humidity	15% - 55% RH
Operating voltage	115/230 V AC, 60/50 Hz, 800W(Approx)
Dimensions (D/W/H)	28" x 39" x 29"
Weight	230 lb/105 Kg

Features and Benefits:

- 30 Slide Open System for creating & optimizing any slide-based protocol
- Compatible with Standard Slides & Micro-Chamber Slides
- Precision Liquid handling for exact Reagent Dispensing
- Cost efficient with minimal Reagent usage
- Green, eco-friendly generates 70% less waste

FISH Processing System: Manual vs NanoVip™

Manual Processing - 33 Complex Steps

NanoVip™ 300 - 3 Simple Steps

[Load](#) → [Click](#) → [View](#)